




THD N-ANO USB

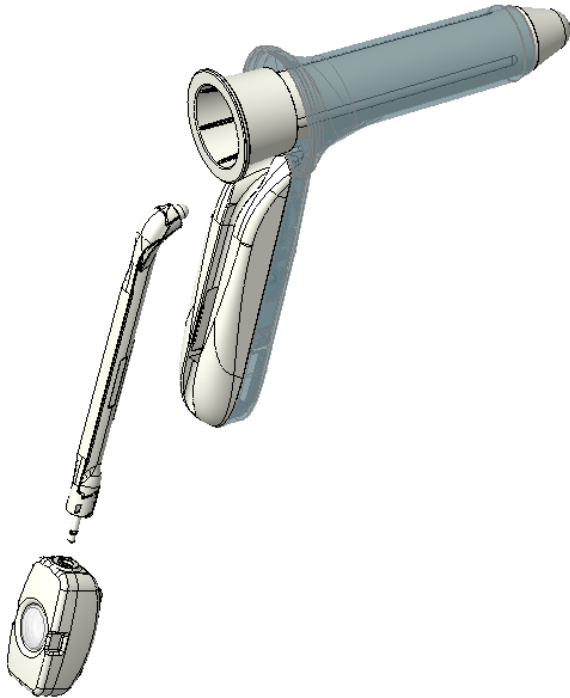

Use

The THD N-ANO USB sales kit includes THD N-ANO BASIC anoscope and its accessories THD PEN-LIGHT and THD LEAD SOURCE USB, intended for the inspection of the anorectal canal.

General Description

THD N-ano BASIC anoscope is a handy and light single use medical device, which presents in the handle a suitable housing useful to supporting THD PEN-LIGHT (THD880070) or THD OPTI-FIBER (THD880012) lighting systems. Available in two dimensions (N-ANO e N-ANO MAXI) each of which can have TRUNCATED or BEVELLED extremity. The THD PEN-LIGHT accessory is a cold light LED source powered by the source dedicated with rechargeable batteries THD LED SOURCE USB. Reusable by sterilization in an autoclave or by other sterilization procedures (see user manual). The THD LED SOURCE USB accessory connects to the THD PEN-LIGHT.

It is equipped with an activation button which can take on two colors according to the state of charge of the device (for the correct interpretation of the status LEDs see user manual REF 370158). The device is rechargeable via standard USB-C cable to be connected to a suitable PC socket. It can not be sterilized in an autoclave so it must be suitably protected during its use.

Characteristics

- In medical grade polycarbonate
- Complies with the biocompatibility requirements established by ISO 10993
- Complies with the electrical safety requirements established by IEC 60601-1
- Complies with electromagnetic compatibility requirements established by IEC 60601-2
- Atraumatic insertion profile
- Extremely compact
- Weight 40 g
- Class I CE device
- Latex Free

Advantages:

- Small size, light and easy to maneuver
- Ready to use
- Minimally invasive, atraumatic



Technical Data Sheet

REF. 800283

REF. 800284

REF. 800285

REF. 800286

EN

THD N-Ano USB						
Code		Q.ty	Trade Name	Description	N° RDM	EMDN
800283 - 50		50	THD N-ANO USB (TR)	N-Ano USB Truncated Anoscope	137358	N.A.
<i>includes</i>	800301	50	THD N-ANO BASIC (TR)	N-Ano Basic Truncated Anoscope	2385382	G02060301
<i>includes</i>	880070	4	THD PEN-LIGHT USB	THD PEN-LIGHT USB	2386595	Z129004
<i>includes</i>	880063	1	THD LED SOURCE USB	THD LED SOURCE USB	2386607	Z129004
800284 - 50		50	THD N-ANO USB (BE)	N-Ano USB Bevelled Anoscope	137347	N.A.
<i>includes</i>	800302	50	THD N-ANO BASIC (BE)	N-Ano Basic Bevelled Anoscope	2385447	G02060301
<i>includes</i>	880070	4	THD PEN-LIGHT USB	THD PEN-LIGHT USB	2386595	Z129004
<i>includes</i>	880063	1	THD LED SOURCE USB	THD LED SOURCE USB	2386607	Z129004
800285 - 50		50	THD N-ANO USB MAXI (TR)	N-Ano USB MAXI Truncated Anoscope	137349	N.A.
<i>includes</i>	800303	50	THD N-ANO BASIC MAXI (TR)	N-Ano Basic MAXI Truncated Anoscope	2385449	G02060301
<i>includes</i>	880070	4	THD PEN-LIGHT USB	THD PEN-LIGHT USB	2386595	Z129004
<i>includes</i>	880063	1	THD LED SOURCE USB	THD LED SOURCE USB	2386607	Z129004
800286 - 50		50	THD N-ANO USB MAXI (BE)	N-Ano USB MAXI bevelled Anoscope	137348	N.A.
<i>includes</i>	800304	50	THD N-ANO BASIC MAXI (BE)	N-Ano Basic MAXI Bevelled Anoscope	2385668	G02060301
<i>includes</i>	880070	4	THD PEN-LIGHT USB	THD PEN-LIGHT USB	2386595	Z129004
<i>includes</i>	880063	1	THD LED SOURCE USB	THD LED SOURCE USB	2386607	Z129004